



Your Safety comes First!


Non-Reinforced SMS

Engineered for demanding situations, providing protection in moderate risk environments and extra comfort when wearing it. It's intended to be used, for example, during long fluid intense procedures, surgery, when pathogen resistance is needed, or infectious diseases are suspected (non-airborne).

Our medical protective apparel has been engineered for all those capeless heroes that fight every day in the front line to save lives. Our main goal is for you and your patients to be safe in every step of the way.



SURGICAL GOWN SPECS

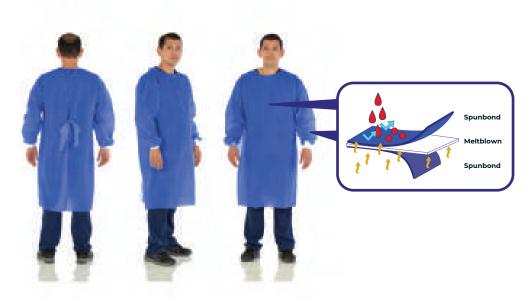

AAMI Level 3 Surgical Gown

510K: K922949 Verasteril-SG3

Protection level: Moderate risk of exposure.

Description: Comfortable, strong, durable and extra refreshing SMS fabric. With nylon stockinette cuffs.

Color: Royal Blue

Size: M, L, XL, XXL, XL-XLONG, XXL-XLONG, XXXL-XLONG

TESTED

Because your safety comes first! We are constantly testing our productos in an **AAMI** approved independent laboratory that certifies our quality compliance.






TESTING THE GOWN

AAMI standards require a Level 3 surgical gown to pass AATCC 42 and AATCC 127 test in the critical zones.

What areas are included in the critical zone?

There are typically 4 different areas that must pass within the critical zones, including two in the front chest and two on the sleeves. If a product passes these tests in all applicable areas, then it qualifies as an AAMI level 3 surgical gown.



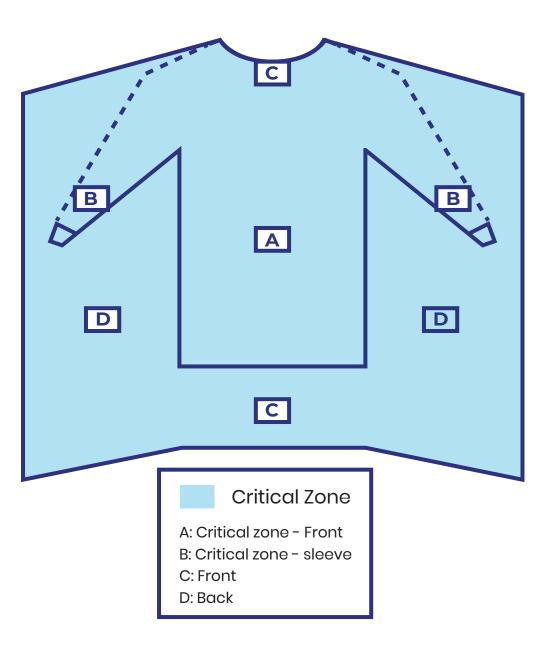

VERASTERIL AAMI LEVEL 3 GOWN				
Test	Standard Used	Liquid Challenge	Result	Expected Barrier Effectiveness
Impact Penetration	AATCC 42	Water	= 1.0 g	Moderate water resistance (resistant to water spray and some resistance to water penetration under constant contact with increasing pressure)
Hydrostatic Pressure	AATCC 127	Water	= 50 cm	



EUROCLYDON LABORATORIES USA LLC



3627 Daintree Park Dr. Katy, TX. 77494-6923

WWW.VERASTERIL.COM
INFO@EUROCLYDONLABS.COM